

Dublitolop[®]



duplicating
material



Dublier-
Gel



gel de
duplication



gelatina
para
duplicar



gelatina da
duplicare

advantages



Dublito[®]

The duplicating material with a difference

The result of the latest Dentaurum research. The reversible high-precision duplicating material is an opaque gel for all phosphate and alcohol based investments which ensures a highly accurate fit when using the CoCr model casting technique and the resin injection technique for dentures as well as for combination work. Dublito[®] can be remelted up to ten times at no loss of quality.

High definition

Optimal contrast due to the new light blue, opaque colour

Extraordinary tensile strength

High stability while unmolding

Excellent fit and accuracy

Best working results

Excellent thermal stability

Cost savings due to long shelf life

No discolouration of the plaster models

Always presentable models

Biodegradable

No disposal problems


Smooth surfaces

Optimum working conditions


Suitable for the microwave oven

Time- and cost saving

Tear resistance chart of Agar Duplicating Materials (original values and after melting 30 times) according to ISO/CD 14356 Dentistry-Duplicating-Material

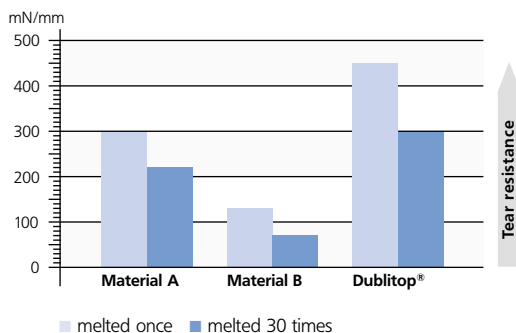

Tensile strength chart of Agar Duplicating Materials (after melting once) according to DIN 53504 testing of rubber and elastomers

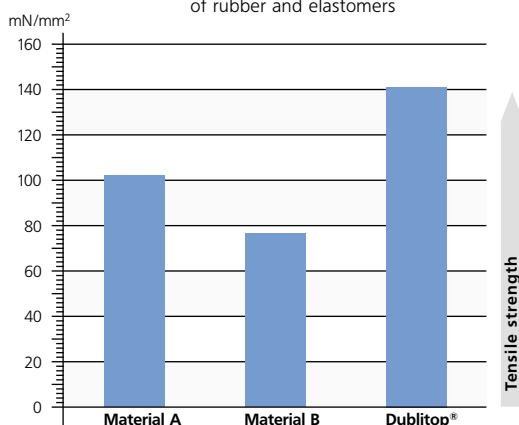

Package unit: Dublito[®]
Order No. 108-610-00
(5kg)



Dublito[®] micro

The innovative portioned package for the microwave oven

Dublito[®] micro, which has the same material characteristics as Dublito[®] (5 kg) has been specially developed for laboratories using relatively small quantities of duplicating material and where an automatic duplicating unit is not economical. Dublito[®] micro now allows all laboratories to benefit from the advantages of Dentaurum high-precision duplicating materials. The portioned pack size makes it easy to re-melt the material as required. A microwave oven is used to prepare the product quickly and in the right quantity without destroying any of the material's characteristics. Dublito[®] micro is outstanding for its dimensional accuracy and its ability to reproduce even the finest details. It can also be stored for long periods of time and is biodegradable. If handled properly, every beaker Dublito[®] micro can be melted at least 10 times.

Gentle melting process preserves material qualities

Incredibly fast application

Minimal technical effort

Easiest handling

Smooth surfaces

Package unit: Dublito[®] micro
Order No. 108-615-00
(4 x 350 g)



Cost effective



Dublito[®] micro Set

With electric mixing unit

An accessory package specially developed to make your work even more convenient and economical. The practical mixer with mixing vessel and digital signal thermometer, is very easy to use and guarantees smooth and even mixing of the duplicating material. This unit allows the molten duplicating material to be cooled down to working temperature in only 10 minutes. When the correct temperature is reached, a signal is displayed. This makes your work easier and saves you time. Quite simply, the Dublito[®] micro Set is yet another highly practical tool for your day-to-day work.



Package unit: Dublito[®] micro Set
Order No. 108-620-00

consists of: Dublito[®] micro (4 x 350 g)
Dublito[®] mixer
Signal thermometer

A good feeling!

D
DENTAURUM

Turnstraße 31 · D-75228 Ispringen · Germany · Telephone +49 72 31/803-0 · Telefax +49 72 31/80 32 95
www.dentaurum.com · E-Mail: info@dentaurum.de