

LUBREX™

Auto Cleaning and Lubricating System



LUBREX™

- Automatic cleaning, purging and lubricating system.
- Ergonomic front-loading, easy access design.
- Can be used for all high speed, low speed, air scaler, and air motor handpieces.
- Up to 4 handpieces can be processed simultaneously.
- Quick release air connector.
- Automatic fan and filter for vapor removal.
- Large reservoir capacity with visible oil level indicator.
- Unit is available in 100V to 240V 50/60 Hz operation.


ISO 9001
ISO 13485



DENTAMERICA®

LUBREX™

Auto Cleaning and Lubricating System

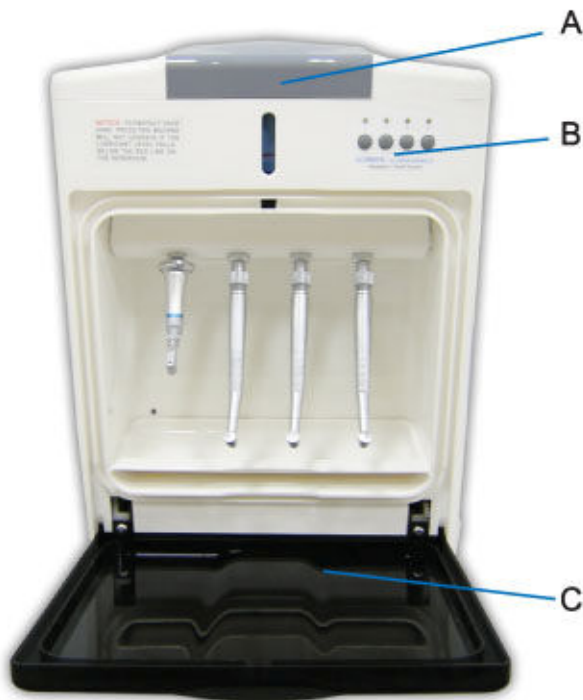

A. Lubricant reservoir



The lubricant reservoir contains Lubrex brand lubricant made from natural mineral oil that is nontoxic, nonflammable, and nonallergenic.

B. Control Panel



The control panel allows for up to 4 handpieces (3 high speed + 1 low speed) to be processed together or individually for maximum efficiency.

C. Front Panel



Innovative design angles the handpiece connectors outward when the front panel is opened.

SPECIFICATIONS

Dimensions	W330 x H430 x D430(mm)
Lubricant tank	300ml
Net Weight	13 kg
Air Pressure	7-10 bar
Main Voltage	100V to 240V
Frequency	50/60HZ
Fuse	1A
Connector types:	
Low-speed handpiece	E-type
High-speed handpiece	Midwest(4-hole) or Borden(2-hole)

DENTAMERICA®

18320 Bedford Circle
P.O. Box 3200
Industry, CA 91744 U.S.A.
Phone: (626) 912-1388
Fax: (626) 913-0510
www.dentamerica.com

DISTRIBUTOR

DENTAMERICA®