

Blendex™

Automatic Alginate Mixer



ISO 9001
ISO 13485



Blendex™

- Automatic trouble-free alginate mixing at the touch of a button.
- Time saving process takes only seconds to complete.
- Consistent, smooth and bubble-free mixing each and every time.
- Features automatic timer with audible tone.
- 3 customizable memory presets available.
- LED display shows mixing time and memory presets.
- Automatic door safety lock ensures safe operation.
- Unit is available in 100/110V 50/60 Hz, or 220/240V 50 Hz operation.

DENTAMERICA®

Blendex™

Automatic Alginate Mixer



① Add alginate and water to the mixing cup.






② Place the mixing cup inside well.





③ Push the START button for a perfectly smooth Blendex™ mix.

3 customizable memory presets available (settings can be changed)

Cups	Mixing Time (sec)
	8
	10
	12

SPECIFICATIONS

Dimensions	W205 x D245 x H295 (mm)
Weight	18 kg
Voltage	AC220V~230V (50/60 Hz) AC100V~120V (50/60 Hz)
Maximum Power	400W
Timer	1-16 sec, adjustable by time
Memory	3 (installed settings: 8, 10, 12 sec)
Frequency	3,600 RPM

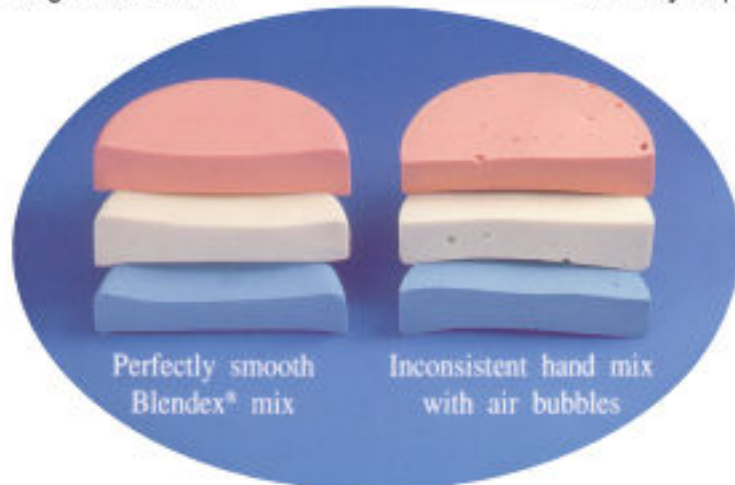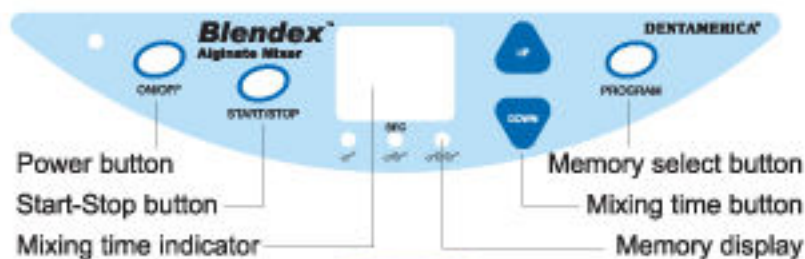

DENTAMERICA®

18320 Bedford Circle
P.O. Box 3200
Industry, CA 91744 U.S.A.
Phone: (626) 912-1388
Fax: (626) 913-0510
www.dentamerica.com

DISTRIBUTOR

DENTAMERICA®