



AIRVAC 2002 –  
**High-performance  
maintenance-free**

The **AIRVAC 2002** is a universal vacuum mixer for use in dental technology. It can be used for mixing any standard investment, plaster or silicone.

A feature of the integrated industrial motor is its very high torque. The high-performance vacuum pump incorporated in the unit is maintenance-free.


## Advantages of the AIRVAC 2002

- very easy to operate
- precise positioning of the mixing shaft contributes to a long service life of the mixing motor
- compact
- automatic regulation of the vacuum provides increased reliability
- the mixing motor remains clean and is easily serviced due to a supporting filter and a fine filter in conjunction with the maintenance-free vacuum pump
- direct connection between mixing bowl and vacuum pump – no disturbing exterior hose
- mechanical coupling locks the mixing bowl in position
- a choice of 3 mixing bowls (900 g, 500 g and 120 g)

## Technical data

### AIRVAC 2002

Mains voltage	230 V
Mains frequency	50 Hz
Average power consumption	250 W
Fuse	T3,15 A
Min. vacuum pressure	100 mbar
Mixer	368 rpm
Dimensions [W x H x D]	150 x 325 x 212 mm
Weight (Wall unit)	11,5 kg

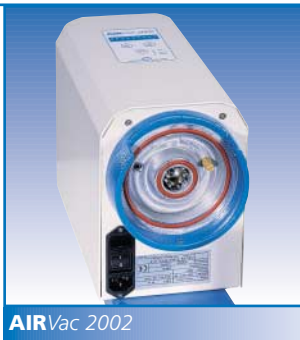

AIRVac 2002



Mixing bowls



Unit base

## Package contents

### AIRVAC 2002

		Order No.
1	AIRVac 2002 vacuum mixer	230 V / 50 Hz 095-072-00
		115 V / 50 Hz 095-073-00
1	Mixing bowl, medium	
1	Instructions for use with EC Declaration of Conformity	
1	Drilling template	
1	Fixing set for mounting on the wall	

## Accessories

### AIRVAC 2002

	Order No.
Unit base	095-075-00
Mixing bowl, large	095-062-00
Mixing bowl, medium	095-063-00
Mixing bowl, small	095-064-00

*A good feeling!*

**D**  
DENTAURUM

Turnstraße 31 · 75228 Ispringen · Germany · Telephone +49 72 31/803-0 · Fax +49 72 31/80 32 95  
www.dentaurum.com · E-Mail: info@dentaurum.de